

Table 1: Comparative features of hand-held and Universal serial bus Dermatoscopes

Feature	Hand-held dermatoscope	Videodermatoscopes
Popular dermatoscopes	DermLite, Heine dermatoscopes	Dinolite Edge Series; Dinolite Basic; Firefly DE300
Technical specifications		
Magnification offered	10x-20x for most devices Heine IC1 offers 12x-40x; Heine NC1 offers 6-9x Additional magnification (upto 50x) possible only through the attached smartphone	10x-220x for most devices Even higher magnifications are available (500-1000x); however, they are not relevant to dermatology These dermatoscopes can visualize the same lesion at two different magnifications (lower magnification at 30-50x and higher magnification at 180-200x)
Power source	USB rechargeable batteries for most devices	No charging required. It comes live when connected to the viewing device
Viewing mode offered	Image can be viewed by the operator only; cannot simultaneously show to the patient Additional screen mirroring applications may be required for the same	Image is viewed on a computer screen, hence, can be simultaneously viewed by patients/ others
Battery life	0.5-8 h depending on the device used	Not applicable
Optical resolution	Very good	Low cost devices have slightly inferior optical resolution; whereas high- end models have a resolution comparable to handheld devices
Polarized and nonpolarized modes	Available in most devices. Toggling possible with a switch	Available in most devices. Toggling possible with the turn of a ring
Contact/noncontact	Most handheld devices are contact devices i.e., the lens comes in direct contact with the lesion being inspected Newer devices are coming up with removable contact face-plates. This helps avoid contact with potentially infected lesions as well as blanching of the vascular structures	These are noncontact devices. Only the rim of the spacer comes in contact with the skin surface. The lens is held back at a distance. It is possible to change the spacer; to sterilize it separately; and also, to cover it to prevent cross-infection
Accessories required		
For viewing	None	Requires a computer/laptop screen to view
For capturing images	Smartphone/camera with/without a connector (Magneticconnect)	Images can be captured on the viewing device
Compatible devices	Smartphones including android phones and iphones; ipads; digital cameras etc	Laptops or computers with a USB port
Versatility of the device		
Pigmentary evaluation	+++	++
Vascular evaluation	++	+++
Trichoscopy	+++	+++
Onychoscopy	++	+++
Nail fold capillaroscopy	+; Blanching of the capillaries is an issue	+++
Entomodermoscopy	++	+++
Mucoscopy	+	+++
Oculoscopy	+	+++
Hidroscopy	+	+++

USB: Universal serial bus

various indications. Table 2 summarizes the common brands available in the market.

For a beginner, especially a post-graduate student, buying an expensive instrument like a dermatoscope is an important investment to make. With limited knowledge and several choices, it is not always easy to make the right decision. This uncertainty prompts many to buy a cheaper dermatoscope in the initial phase of their learning. However, poor resolution and picture quality fail to do justice to dermatoscopy, and may

dampen one's interest. Therefore, it is advisable to research available instruments and invest in the best one possible.

The best dermatoscope to buy is the one which best suits the user's personal needs in terms of the type of skin lesion to be visualized commonly. This is further based on the type of practice, time available and the ease of carrying accessory equipment for image capture and storage [Figures 1 and 2]. Hand-held dermatoscopes offer low magnification but excellent optical resolution and images, but most devices

Table 2: Commonly available dermatoscope brands

Device	Company	Imaging technology
Hand-held dermatoscopes		
DermLite Cam	3Gen	Built-in digital camera
DermLite DL1	3Gen	Mobile device
DermLite DL200, DermLite DL3/DL4	3Gen	Mobile device or Digital Camera
DermLite FOTO System, DermLiteFoto II Pro	3Gen	DSLR
VEOS HD1/HD2	Canfield	Mobile Device
VEOS SLR	Canfield	DSLR
NC2 Dermatoscope	Heine	Mobile Device
DELTA 20T	Heine	DSLR
iC1	Heine	Mobile Device
Illuco IDS-1100	Dermoscan	Mobile Device or digital camera
Video dermatoscopes		
Medicam 1000	FotoFinder	Built-in video dermoscopy
DermoGenious Ultra	Dermoscan	Built-in video dermoscopy
MoleMax	Derma Medical Systems	Built-in video dermoscopy
Dino-Lite Edge Series, Dino-Lite Premier Series	AnMo Electronics Corporation	Built-in video dermoscopy
FireFly DE300, FireFly DE350	FireFly	Built-in video dermoscopy
Hi-Scope	Hirox	Built-in video dermoscopy
VideoCap	D.S. Medica	Built-in video dermoscopy
Horus	Adamos.r.l	Built-in video dermoscopy
Optilia	Optilia instruments	Built-in video dermoscopy
VivaCam	Caliber Imaging and Diagnostics	Built-in video dermoscopy
C-Cube	Pixience	Built-in video dermoscopy

need to be in contact with the area to be examined, which may limit their use, especially in the current pandemic era. Other advantages of hand-held dermatoscope include ease of use and easy connectivity to the smartphones. For daily practical needs, a hand-held scope might be as good as video one or even better given the practicality of use, while for research purposes, a videodermoscope might be better. Videodermoscopes are more versatile devices which allow visualization of odd sites like nail, eye and mucosae with ease.⁶⁻⁹ Features like magnification, resolution, ultraviolet light, pigment boost, measurement and calibration add to the versatility of these dermatoscopes. The availability of

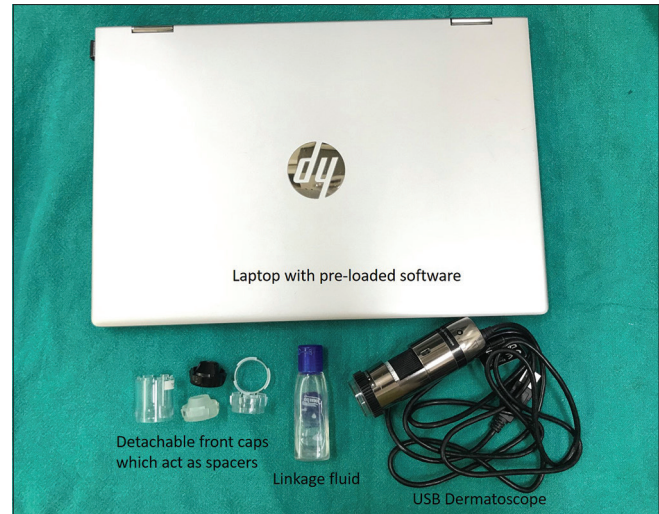


Figure 2: A universal serial bus dermatoscope with accompanying requirements for capturing and storing images

cross-polarization should always be ensured in the device bought. It is recommended to hold and use the device before investing a major amount, so that individual comfort and ergonomics with the procedure are ensured. Whatever type of dermatoscope one starts off with, there always exists the possibility of switching to another type to discern more features or to capture better images.

Financial support and sponsorship
Nil.

Conflicts of interest
There are no conflicts of interest.

References

1. Kaliyadan F. The scope of the dermatoscope. *Indian Dermatol Online J* 2016;7:359-63.
2. Micali G, Lacarrubba F. Dermatoscopy: Instrumental update. *Dermatol Clin* 2018;36:345-8.
3. Jakhar D, Grover C, Kaur I, Sharma S. Response to “Dermatoscopic features of lichen nitidus”. *Pediatr Dermatol* 2019;36:179.
4. Lacarrubba F, Dall’Oglio F, Rita Nasca M, Micali G. Videodermatoscopy enhances diagnostic capability in some forms of hair loss. *Am J Clin Dermatol* 2004;5:205-8.
5. Steeb T, Wessely A, Niesert AC, Ruzicka T, von Braunmühl T, Berking C, *et al.* Patient attitude towards videodermatoscopy for the detection of skin cancer: A cross-sectional study. *Oncol Res Treat* 2019;42:319-25.
6. Jakhar D, Grover C, Singal A. Nailfold capillaroscopy with USB dermatoscope: A cross-sectional study in healthy adults. *Indian J Dermatol Venereol Leprol* 2020;86:33-8.
7. Jakhar D, Grover C. Innovative modification of the USB dermatoscope for mucoscopy. *J Am Acad Dermatol* 2018;78:e3-e4.
8. Jakhar D, Grover C. Universal serial bus dermatoscope as an oculoscopy tool. *J Am Acad Dermatol* 2018;78:e139-e140.
9. Jakhar D, Kaur I. Hidroscopy: *In vivo* videodermoscopy of the sweat glands. *Skin Res Technol* 2019;25:410-1.